

SEGUFIX®-Shoulder Restraint



Item 2203 / Item 2303 / Item 2403

Consisting of:

- 1 SEGUFIX®- Shoulder Restraint
- 2 SEGUFIX®- Locks (keys not included)

Sizes: XS / S / M / L / XL

Versions: basic version / r

Component part: item 2103
(locks and keys not included)



The SEGUFIX®-Shoulder Restraint (red label) is part of the SEGUFIX®-Bed Restraint System and made from a cotton/rayon mix.

Note: SEGUFIX®-Standard required (see image)

3 Versions → SEGUFIX®-Locks:		
1206 / 1208	1306 / 1308	1406 / 1408
Item 2203	Item 2303	Item 2403

Table of Contents	Page
SEGUFIX®-Shoulder Restraint – Instructions	2 - 4
Overview table: all versions, sizes, measurements and weights	5
SEGUFIX®-Locks and Keys	6 - 7
SEGUFIX®-Safety Information	8 - 11
Instructions for Care	12
Imprint	12

Versions – Item no. / identification	
Item no. / Version	Product features
...	Basic version
.../r	Eyelet strap reinforced for higher demands


SEGUFIX®-Shoulder Restraint

3 Locking Systems:

Item 2203 / 2303 / 2403

Intended use:

- Prevents the patient from sliding up out of the SEGUFIX®-Standard
- Physically restraining patients (e.g. pathology or disability that requires restraint in a defined position in bed to restrict patients from harming themselves or care givers, e.g. anxiety disorders)

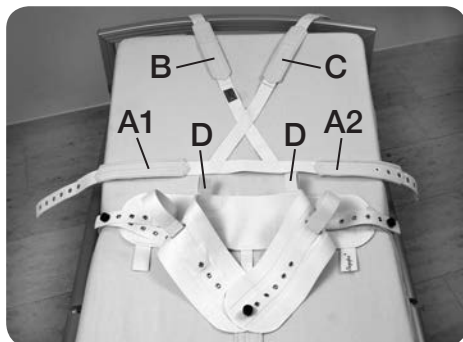

Instructions for use item 2203 / 2303 / 2403 and 2103



Locks and keys see pages 6 and 7

First install the SEGUFIX®-Standard to the bed
(please refer to instructions for SEGUFIX®-Standard)

How to attach the SEGUFIX®-Shoulder Restraint to the waist belt:

The SEGUFIX®-Shoulder Restraint consists of three connected eyelet straps:
The underarm strap (A1/A2) and the shoulder straps (B and C).

- Lay the underarm strap (A1/A2) parallel to the waist belt of the Standard onto the bed. The metal loops (4) on the shoulder straps (B and C) should be facing the mattress.
- Pull the ends of the underarm strap (A1/A2) through the respective loops (D), located on the waist belt of the SEGUFIX®-Standard, towards the edge of the bed.

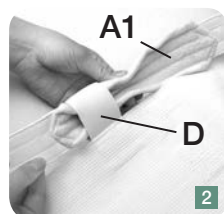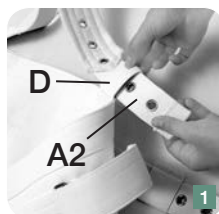



- The Shoulder Restraint is now attached to the waist belt of the SEGUFIX®-Standard and ready to be fastened to the patient.

How to fasten the SEGUFIX®-Shoulder Restraint to the patient:

- Insert the metal pin (1) through the respective eyelet of the pin flap on the shoulder strap (C).



When using the SEGUFIX®-Standard with Thigh Straps:

- Pull both shoulder straps (B and C), that

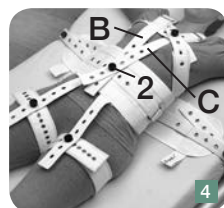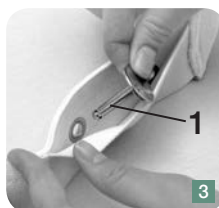


are pointed towards the head of the bed, tightly over the shoulders of the patient and lay them over the pin of the open lock (2) located in the middle of the waist belt of the Standard.

- Close the lock (2) (fig. 4).
The metal loops (4) on the shoulder straps are now facing up.
- Pull the end of the underarm strap (A1) under the right armpit of the patient, insert it through one of the three slots (E) on the shoulder strap (B), and lay the respective eyelet onto the pin of the lock (1), so that the underarm strap (A1) fits tightly.
- Next, pull the other end of the underarm strap (A2) under the left armpit of the patient, and lay it likewise tightly fitting onto the pin of the lock (1).
- Close the lock (1).
- Pull on each lock to test for a secure hold. Then lock bed side rails in an up position.

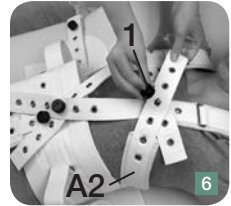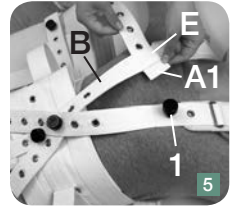

When using the SEGUFIX®-Standard with Crotch Strap (all versions) or the SEGUFIX®-Standard with Thigh Straps (versions “r” or “f”):

- Insert the pin of another lock (3) (included in the delivery) through one eyelet of the crotch strap or the waist belt, pull both shoulder straps (B and C) tightly over the shoulders of the patient and lay them over the pin of the additional lock (3).
The metal loops (4) on the shoulder straps are now facing up.
- Close the lock (3).
- Pull the end of the underarm strap (A1) under the right armpit of the patient, insert it through one of the three slots (E) on the shoulder strap (B), and lay the respective eyelet onto the pin of the lock (1), so that the underarm strap (A1) fits tightly.
- Next, pull the end of the underarm strap (A2) under the left armpit of the patient, and lay it likewise tightly fitting onto the pin of the lock (1).
- Close the lock (1).
- Pull on each lock to test for a secure hold.
Then lock bed side rails in an up position.

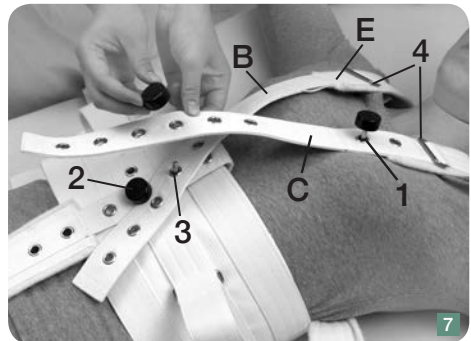



Important:

The underarm strap (A1/A2) must be secured just under the armpit, otherwise the patient could escape. Make sure that none of the straps hinder or restrict the patient's neck or throat. The straps must be tight but should not restrict breathing.

How to release the restraint and remove it:

See page 7 on how to open the respective locking systems.

- Open the locks (1 and 3), pull straps off of the pins and remove them from the loop (E).
- Then remove the Shoulder Restraint from the loops (D) of the waist belt of the Standard, and remove it from the bed.

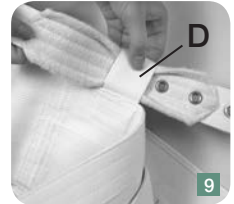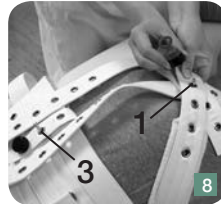

Follow SEGUFIX®-Safety Information and Warnings on pages 8 - 11

WARNING: Read the Safety Information prior to using the SEGUFIX®-Shoulder Restraint.

- The SEGUFIX®-Shoulder Restraint can only be used in connection with the SEGUFIX®-Standard – all versions with locking systems (no Velcro® products).
- Locks under tension may be difficult to open, or not open at all. Therefore use a maximum of 4 belt layers per lock.

Version “r” or “f” products are thicker. These products allow for a maximum of 3 belt layers per lock. When the SEGUFIX®-Shoulder Restraint is used in connection with the SEGUFIX®-Standard with Crotch Strap (all versions) or Standard with Thigh Straps (versions “r” or “f”), use an additional lock (included with the delivery).

When using the Patent Twist Lock (item no. 1306) use no more than 4 belt layers.

When using the Push Lock (item no. 1406) use no more than 3 belt layers.

- The SEGUFIX®-Shoulder Restraint does not prevent the patient from sitting-up in bed. To prevent the patient from sitting-up in bed, the SEGUFIX®-Shoulder Additional Belt must be used.
- Pacemaker warning: Ensure a safety distance of 10 cm (4 inches) between the SEGUFIX®-Magnetic locking system / the SEGUFIX®-Patent Twist locking system and the pacemaker (alternatively use the Push Lock).
- Note that attachment straps and moving mechanics of the bed or mattress may interfere with each other, causing damage, loss of function and injury.
- For more product info visit us online at www.segufix.com or www.segufix-shop.com




Example of use:
SEGUFIX®-Shoulder Restraint and
SEGUFIX®-Shoulder Additional Belt,
prevents the patient from sitting up in
bed (ensure continuous monitoring)

SEGUFIX®-Shoulder Restraint

Size	Item	kg	Item	kg	Item	kg	Component part	kg
XS	2203	0,24	2303	0,28	2403	0,28	2103	0,22
S	2203	0,48	2303	0,52	2403	0,52	2103	0,44
M	2203	0,52	2303	0,56	2403	0,56	2103	0,50
L	2203	0,56	2303	0,60	2403	0,60	2103	0,54
XS	2203/r	0,26	2303/r	0,30	2403/r	0,30	2103/r	0,24
S	2203/r	0,47	2303/r	0,51	2403/r	0,51	2103/r	0,46
M	2203/r	0,55	2303/r	0,59	2403/r	0,59	2103/r	0,53
L	2203/r	0,59	2303/r	0,63	2403/r	0,63	2103/r	0,57
XL	2203/r	0,63	2303/r	0,67	2403/r	0,69	2103/r	0,61

Size	Width of strap		Length of strap			
	mm	inches	Strap a mm	Strap b mm	Strap a inches	Strap b inches
XS	35	1.4	600	770	23.6	30.3
S	45	1.8	860	1190	33.9	46.9
M	45	1.8	1010	1450	39.8	57.1
L	45	1.8	1170	1650	46.1	65.0
XL	45	1.8	1420	1870	55.9	73.6

Table: approximate measurements – about 8 % shrinkage possible (products are not pre-washed)

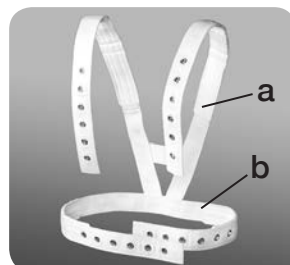

Length of strap: strap a and strap b, see table



SEGUFIX®-Shoulder Restraint
Back view size XS / S



SEGUFIX®-Shoulder Restraint
Back view size M / L / XL



SEGUFIX®-Magnetic locking system black

Item	Consisting of:
1206	1 Magnetic Button black / segufix.com 1 Metal Pin
1208	1 Magnetic Key red
1204	1 Magnetic Button black / segufix.com
1205	1 Metal Pin
1209	1 Magnetic Key red with key ring



Push locking system

Item	Consisting of:
1406	1 Push Button 1 Push Pin
1408	1 Push Lock Key
1404	1 Push Button
1405	1 Push Pin



SEGUFIX®-Patent Twist locking system green

Item	Consisting of:
1306	1 Patent Twist Button green 1 Square Pin
1308	1 Patent Twist Key yellow
1304	1 Patent Twist Button green
1305	1 Square Pin

Intended use of the Lock:

Closing of the SEGUFIX® systems.

Intended use of the Key:

Quick release of the locks.

Caution: Each locking system only works with its own component parts and is not interchangeable with other locking systems.



WARNING: Read the Safety Information prior to using the SEGUFIX®-Locks and Keys.

- Do not submerge locks in any type of liquid. This can cause malfunction.
- Locks and keys must not be laundered.
- Test locks for proper function before each use.
- Examine the head of the pin for damage (wear and tear, rounding off or deformation).
- When closing the lock make sure that the button locks securely onto the pin. You should not be able to pull the button off of the pin.
- When opening locks, ensure movement between the upper and lower points of contact (item 1206). Therefore each lock can hold up to a maximum of 4 belt layers (exception: products with "r" or "f" in the item number have thicker belts. These allow for a maximum of 3 belt layers!). When using the Patent Twist Lock (item no. 1306) use no more then 4 belt layers. When using the Push Lock (item no. 1406) use no more then 3 belt layers.
- If damaged locks fail to open, straps can be cut with scissors or belt cutters.
- For more product info visit us online at www.segufix.com or www.segufix-shop.com

Instructions for use: SEGUFIX®-Locks and Keys

Item 1206 / 1208 (SEGUFIX®-Magnetic locking system black)

Item 1306 / 1308 (SEGUFIX®-Patent Twist locking system green)

Item 1406 / 1408 (Push locking system)

How to open the Locks:

Item 1206 / 1208

- Set the red Magnetic Key with a light push on top of the black Magnetic Lock.
- Lift both the Magnetic Key and the Magnetic Button off of the Metal Pin (fig. 1).



Item 1306 / 1308

- Set the yellow Magnetic Twist Key on top of the green Magnetic Twist Lock.
- Unlock the Patent Twist Lock by turning the Patent Twist Key, and then lift both the Patent Twist Key and the Patent Twist Button off of the Square Pin (fig. 3).



Item 1406 / 1408

- Push the two pins of the Push Lock Key into the holes on top of the Push Button to unlock.
- Lift both the Push Lock Key and the Push Button off of the Push Pin (fig. 2).



How to close the Locks:

- Set the respective button with the opening onto the pin so that it locks in place.
- Now test the lock for a secure hold: You should not be able to pull the button off of the pin.

Caution: Each locking system only works with its own component parts and is not interchangeable with other locking systems.

Follow SEGUFIX®-Safety Information and Warnings on pages 8 - 11

The following instructions for use are only applicable to original SEGUFIX® products. Do not use the SEGUFIX®-Restraint and Positioning System until you have read and understood the instructions for use.

Only use up-to-date instructions and keep them in a safe and readily available location. Product and instruction must match. Failure to use this product properly can cause serious injury or death. Your patient's safety depends on your using the SEGUFIX® system correctly.

1. Requirements for use

It is the user's responsibility to comply with all national laws. Always consider restraining protocols of your work place and current legislation. The instructions for use must be adhered to.

SEGUFIX® products may only be used by individuals with proper training and knowledge of the product. Legal provisions and doctor's orders must be strictly observed! The attending physician must ensure, that the patient's medical condition allows the application of the SEGUFIX® system.

2. Intended use

SEGUFIX® products may only be used for the medical purposes, which they were designed for.

3. Remove dangerous objects

Remove all objects from the reach of the patient, which could lead to injury or may damage the belts (e.g. glasses, jewelry and sharp objects).

4. Inspect products before use



Inspect before each use for broken stitches, torn, worn out or frayed parts and damaged eyelets. Damaged products could malfunction and cause injury. Do not use damaged or worn out products.

5. When opening, locks may not be under tension

When opening locks, ensure movement between the upper and lower points of contact (item 1206). Therefore each lock can hold up to a maximum of 4 belt layers (exception: products with "r" or "f" in the item number have thicker belts. These allow for a maximum of 3 belt layers!). When using the Patent Twist Lock (item 1306) use no more than 4 belt layers. When using the Push Lock (item 1406) use no more than 3 belt layers.



6. Always attach Side Positioning Straps

The SEGUFIX®-Standard must always be used with the Side Positioning Straps attached. Side Positioning Straps prevent the patient from helicoptering in bed or from getting over the edge of the bed causing compression of the chest, possibly resulting in death.

In some cases it is possible that patients get over the edge of the bed even with the Side Positioning Straps attached and the side rails in the up position. In these cases additional restraint measures must be taken.

Important: Insert the Side Positioning Straps through the respective loops (according to instructions) before securing them to the bed belt.

7. Beds, chairs and stretchers


- Only use SEGUFIX® products on beds, stable chairs and stretchers that allow the product to be attached to, in accordance with the instructions.
- Products must be attached in such manner, that they cannot get out of place.
- Note that straps and moving mechanics of beds, stretchers or chairs (or alternating pressure pads and the like), may interfere with each other causing damage, loss of function and injury.
- The attachment points must be free of sharp edges.
- SEGUFIX® products do not prevent beds, chairs or stretchers from tipping over or from unintended movement.

Helpful hints:

European standards for beds: DIN EN 60601-2-52

European standards for wheelchairs: DIN EN 12182 and DIN EN 12183

8. All side rails must be in the up position

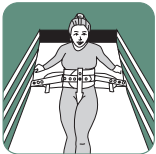

To prevent accidents always lock side rails in the up position (exception: 5-point-restraint).

Caution: Do not use with parted side rails except with a gap protector that effectively prevents the patient from slipping through the gap. An open gap could allow the patient to slip through even with the restraint fastened possibly causing compression of the chest.

9. Monitor patients in restraints



Patients in restraints require frequent monitoring.

Aggressive, agitated or restless patients and patients in danger of aspiration require constant monitoring.

10. Test closure-systems prior to each use



Check SEGUFIX®-Magnetic Locks black/segufix.com, SEGUFIX®-Patent Twist Locks green, Push Locks, Snap Locks, Velcro® Fasteners and Buckle Locks for a secure hold.

Do not submerge locks in any type of liquid. This may cause malfunction. If damaged locks fail to open, straps can be cut with scissors or belt cutters. Check the head of the pin for signs of wear and tear (rounding off

or deformation).

11. Pacemaker Warning



A distance of 10 cm (4 inches) must be maintained between a patient's pacemaker and magnetic locks and keys to prevent the pacemaker from possible malfunction (alternatively use Push Locks).

For possible interference of strong magnets with other patient implants please refer to manufacturer's instructions of these products.

12. Test products for correct attachment

Frequently check products for correct attachment. The straps must be secured to the frame of the bed (or the folding part of the frame of the bed) in such a way, that they cannot slide out of place.

13. Only use original SEGUFIX® products

For safety reasons SEGUFIX® products must not be altered or be combined with third party products.

14. Use of SEGUFIX® products in vehicles



The SEGUFIX®-Transport Belt System does not replace the vehicle's legally required seat belt system. SEGUFIX®-Transport Belts may only be used as additional support or for restraint of a patient. Make sure that the patient can be rescued quickly in case of emergency (e.g. clearly visible and handy belt cutters).

15. Secure but not too tight

The restraints must be fastened tightly, but not interfere with breathing or blood circulation. Loosely fastened products may jeopardize the patient's safety. Frequently check products for correct position.

16. Storage

To prolong product life store belts and locking systems in dry and dark conditions at about 20 °C.



17. Flammability of SEGUFIX®-Bed Restraints



According to DIN EN 71-2 section 4.3 non-flammable.

According to DIN 75200 the burning rate is zero.

However, products may not come in contact with cigarettes or be exposed to open flame, as this could permanently damage the products, and this could lead to injury.

18. Correct sizing

Choosing the right size is of essential importance. Make sure that the size matches body and waist sizes (refer to sizing tables). Products that are too small or too large compromise the patient's comfort and safety.

19. Disposal

Carton and plastic wrapping can be disposed of in recyclable disposal, defective products in household garbage.

20. Declaration of conformity

For the declaration of conformity please visit:

https://www.segufix.com/conformity_en.php

21. Expert recommendation to prevent accidents

(Prof. Dr.-Ing. U. Boenick, TU-Berlin)

1. Always use Side Positioning Straps
2. Lock non-parted side rails in the up position
(exception: five-point-restraint)
3. The use of additional restraint elements may be necessary to eliminate bed egress.

22. Ensure prior to each use

- Test locks for a secure hold
- Visual inspection of belts, straps and eyelets
- Ensure sufficient medical indication
- Check legal requirements

Care Instructions

- **SEGUFIX®-Restraints** (white) are made of a cotton/rayon mix and are washable in hot wash cycles up to 95 °C. Lower temperature and laundry nets will prolong product life. Avoid methods that might bend metal parts of the system (e.g. laundry line, dry pressing, rotary irons), eyelets may be damaged.
Not pre-washed, up to 8 % shrinkage. Store in a dry place.
Detergents: non-aggressive, no fabric softener, no bleach.
Dryer: delicate cycle at low temperature and preferably in the laundry net.
Disinfection: may be disinfected using chemothermal disinfection washing procedure.
Overdosing can damage the product.
- **SEGUFIX®-Transport Belts** (coloured/black) are made of polyester and are washable in warm wash cycles up to 30 °C.
Avoid methods that might bend metal parts of the system (e.g. laundry line, dry pressing, rotary irons), eyelets may be damaged. Store in a dry place.
Detergents: non-aggressive, no fabric softener, no bleach.
Dryer: delicate cycle at low temperature and preferably in the laundry net.
- **SEGUFIX®-Locks and metal pins**: remove all locks and metal pins before washing or drying. Laundering will damage locks and possibly cause them to fail. Locks may be cleaned with a damp cloth. Never submerge locks in any type of liquid.
- **Velcro® Fastener**: fasten before washing and drying to prevent lint buildup. Occasionally use a stiff brush to remove lint from Velcro® Fastener. Even with normal use, Velcro® material has the tendency to collect lint. This can reduce the grip strength. Before each use, test grip strength of the Velcro® Fastener.
- **Durability**: depends on the intensity of use and is prolonged by gentle laundry procedures (detergent/laundry net).

Edition 15 · 2019-06

© by SEGUFIX®-Bandagen, 2019

The general business and delivery terms apply:

https://www.segufix.com/PDFs/AGBs/Terms_and_Conditions.pdf

Product specifications subject to change.
Deviations from pictures in colour and shape, errors and amendments excepted.

All rights reserved. Reprints or distribution using electronic systems/media, in its entirety or in parts, with prior explicit authorization in writing only.

All requests concerning this must be directed at:

SEGUFIX®-Bandagen

Das Humane System GmbH & Co. KG

Allerbeeksring 33 · 21266 Jesteburg
Germany

Phone: +49 (0) 41 83 500-0

Fax: +49 (0) 41 83 500-200

E-mail: export@segufix.de

www.segufix.com

www.segufix-shop.com