

# SEGUFIX®-Standard with Crotch Strap



The SEGUFIX®-Standard with Crotch Strap (yellow label) is part of the SEGUFIX®-Bed Restraint System and is made from a cotton/rayon mix.

## Item 2221 / Item 2321 / Item 2421

- Consisting of:**  
 1 SEGUFIX®-Standard with Crotch Strap and integrated Side Positioning Straps  
 3 SEGUFIX®-Locks  
 1 SEGUFIX®-Key

**Sizes:** S / M / L / XL / XXL\*  
 (\*size XXL on request)

**Versions:** basic version / b / e / f / r

**Component part:** item 2121  
 (locks and keys not included)

3 Versions → SEGUFIX®-Locks:		
1206 / 1208	1306 / 1308	1406 / 1408
Item 2221	Item 2321	Item 2421

Table of Contents	Page
SEGUFIX®-Standard with Crotch Strap – Instructions	2 - 4
Overview table: all versions, sizes, measurements and weights	5
SEGUFIX®-Locks and Keys	6 - 7
SEGUFIX®-Safety Information	8 - 11
Instructions for Care	12
Imprint	12

Versions – Item no. / identification			
Item no. / Version	Product features	Span width*	
		mm	inches
...	Attachment strap: medium	1200 - 1450	47.2 - 57.1
... /b	Attachment strap: short	1000 - 1200	39.4 - 47.2
... /e	Attachment strap: long	1450 - 1700	57.1 - 66.9
... /f	Fleece, for sensitive skin, 100 % pure cotton		
... /r	Eyelet strap reinforced for higher demands		


\*Span width: distance between the attachment points on the bed  
  
 approx. measurements – about 8 % shrinkage possible (products are not pre-washed)

# SEGUFIX®-Standard with Crotch Strap

3 Locking Systems:  
Item 2221 / 2321 / 2421

## Intended use:

- Patient restraint and positioning system in bed
- Physically restraining patients (e.g. pathology or disability that requires restraint in a defined position in bed to restrict patients from harming themselves or care givers, e.g. anxiety disorders)


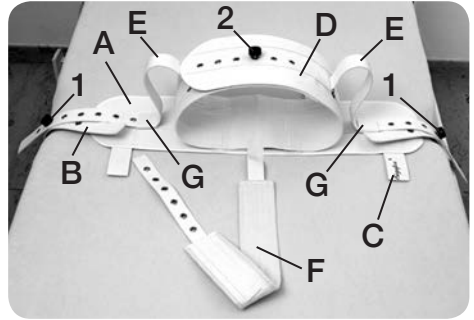


## Instructions for use item 2221 / 2321 / 2421 and 2121

Locks and keys see pages 6 and 7



### How to attach the SEGUFIX®-Standard with Crotch Strap to the bed:

- Insert the pin of the lock (1) through one of the eyelets of the pin flap.
- Lay the wide bed belt (A) with the attachment straps (B) on the mattress at the height of the patient's waist so that the Crotch Strap (F) and the loops (C) are pointed towards the foot of the bed.
- Direct the attachment strap (B) from inside out around the folding part of the bed, without folding part around the bed frame. Install the attachment strap in such manner that it cannot slide along the bed frame. Be sure to pull the attachment strap firmly down and then up around the frame. Then lay the attachment strap tightly with one eyelet over the pin of the lock (1).
- Close the lock (1).
- Install the second attachment strap in the same manner.
- Ensure that the attachment straps are installed very tightly to the bed.

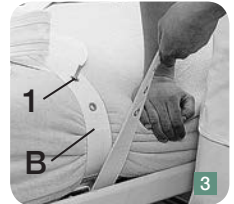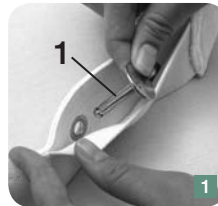

### How to fasten the waist belt to the patient:



- Insert the pin of the lock (2) through one of the eyelets of the pin flap. The metal pin must be in the center of the closed waist belt.
- Wrap the waist belt (D) around the patient's waist, with an eyelet over the pin of the lock (2).

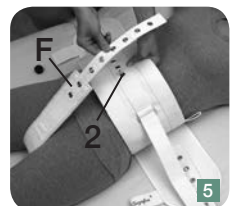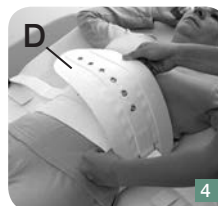

- Pull the Crotch Strap (F) up between the patient's legs towards the waist belt.
- Lay the Crotch Strap (F) with **only one eyelet** onto the pin of the lock (2).
- Close the lock (2).
- The Crotch Strap must fit snugly, but not cause pressure marks.
- The waist belt must fit snugly, but not restrict breathing.  
An open hand should fit between patient and waist belt.

## How to fasten the Side Positioning Straps:


- Insert the Side Positioning Strap (E) through the loop (G).
- Open the lock (1).
- Lay the Side Positioning Strap (E) with an eyelet over the pin of the lock (1).
- Close the lock (1).
- Fasten the second Side Positioning Strap in the same manner.
- The lateral movement can be adjusted by selecting one of the many eyelets according to indication.
- Pull on each lock to test for a secure hold. Then lock bed side rails in an up position.

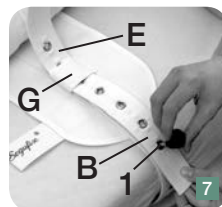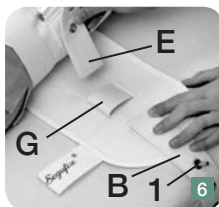

## How to release the restraint and remove it:

See page 7 on how to open the respective locking systems.

- To release the restraint open the lock (2) and lift all belt parts off of the pin.
- Then open both locks (1) of the attachment straps and remove the restraint from the bed.

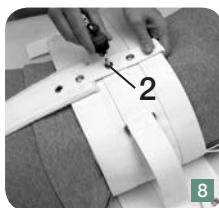

Follow **SEGUFIX®-Safety Information and Warnings** on pages 8 - 11

## Note:

- The SEGUFIX®-Standard with Crotch Strap and Side Positioning Straps does not prevent the patient from sitting up in bed. For this purpose use the SEGUFIX®-Shoulder Restraint in connection with the SEGUFIX®-Shoulder Additional Belt.
- To prevent the patient from sliding up out of the SEGUFIX®-Standard with Crotch Strap, the SEGUFIX®-Shoulder Restraint must be used.
- If “r” and “f” version products are used, each lock can hold a maximum of three belt layers to ensure easy opening of the lock. More than three belt layers put the lock under tension, making it difficult to open or not open at all. When using more than three belt layers, an additional lock is required.
- Loop C: This may be used in connection with the SEGUFIX®-Hand Restraint, see instructions.

**WARNING:** Read the Safety Information prior to using the SEGUFIX®-Standard with Crotch Strap. Improper use can result in serious injury or compression of the chest possibly causing death. Therefore, the following must always be adhered to:

- The SEGUFIX®-Standard with Crotch Strap must never be used without Side Positioning Straps. Otherwise the patient could helicopter in bed or turn over the edge of the bed, possibly resulting in compression of the chest.
- The Site Positioning Straps must always be inserted through the loop (G) of the wide bed belt (fig. 6 and 7, page 3). Exception: five-point-restraint.
- In some cases it is possible that patients get out of bed or over the edge of the bed even with the Site Positioning Straps attached and the side rails in the up position. In these cases additional restraint measures must be taken, e.g. SEGUFIX®-Shoulder-Restraint, three- or five-point-restraint etc.
- Bed side rails must always be locked in the up position. Use only non-parted side rails, or parted side rails with a gap protector (except five-point-restraint).
- The attachment straps (B) of the wide bed belt (A) must be installed in such manner that the wide bed belt lies tightly on the mattress.
- Note that attachment straps and moving mechanics of the bed or mattress may interfere with each other, causing damage, loss of function and injury.
- For more product info visit us online at [www.segufix.com](http://www.segufix.com) or [www.segufix-shop.com](http://www.segufix-shop.com)

## Follow SEGUFIX®-Safety Information and Warnings on pages 8 - 11



Example of use with SEGUFIX®-Magnetic locking system black

Size guide			
Size	Waist girth	Waist belt ①	Bed belt ②
S	520 - 820 mm 20.5 - 32.3 inches	140 mm 5.5 inches	180 mm 7.1 inches
M	800 - 1200 mm 31.5 - 47.2 inches	180 mm 7.1 inches	180 mm 7.1 inches
L	1160 - 1560 mm 45.7 - 61.4 inches	180 mm 7.1 inches	180 mm 7.1 inches
XL	1360 - 1800 mm 53.5 - 70.9 inches	180 mm 7.1 inches	180 mm 7.1 inches
XXL*	1560 - 1980 mm 61.4 - 78.0 inches	180 mm 7.1 inches	180 mm 7.1 inches


\*size XXL on request

Tables: approximate measurements - about 8 % shrinkage possible (products are not pre-washed)

SEGUFIX®-Standard with Crotch Strap									Span width*	
Size	Item	kg	Item	kg	Item	kg	Component part	kg	mm	inches
S	2221	1,06	2321	1,15	2421	1,14	2121	0,98	1200-1450	47.2 - 57.1
M	2221	1,28	2321	1,37	2421	1,36	2121	1,20	1200-1450	47.2 - 57.1
L	2221	1,44	2321	1,51	2421	1,50	2121	1,34	1200-1450	47.2 - 57.1
S	2221/b	1,01	2321/b	1,10	2421/b	1,09	2121/b	0,93	1000-1200	39.4 - 47.2
M	2221/b	1,23	2321/b	1,32	2421/b	1,31	2121/b	1,15	1000-1200	39.4 - 47.2
L	2221/b	1,39	2321/b	1,48	2421/b	1,47	2121/b	1,29	1000-1200	39.4 - 47.2
S	2221/e	1,12	2321/e	1,21	2421/e	1,20	2121/e	1,04	1450-1700	57.1 - 66.9
M	2221/e	1,34	2321/e	1,41	2421/e	1,40	2121/e	1,24	1450-1700	57.1 - 66.9
L	2221/e	1,48	2321/e	1,59	2421/e	1,58	2121/e	1,42	1450-1700	57.1 - 66.9
S	2221/r	1,16	2321/r	1,23	2421/r	1,22	2121/r	1,06	1200-1450	47.2 - 57.1
M	2221/r	1,40	2321/r	1,41	2421/r	1,40	2121/r	1,30	1200-1450	47.2 - 57.1
L	2221/r	1,54	2321/r	1,57	2421/r	1,56	2121/r	1,44	1200-1450	47.2 - 57.1
XL	2221/r	1,66	2321/r	1,71	2421/r	1,70	2121/r	1,56	1200-1450	47.2 - 57.1
S	2221/rb	1,10	2321/rb	1,16	2421/rb	1,16	2121/rb	1,00	1000-1200	39.4 - 47.2
M	2221/rb	1,34	2321/rb	1,40	2421/rb	1,40	2121/rb	1,24	1000-1200	39.4 - 47.2
L	2221/rb	1,48	2321/rb	1,54	2421/rb	1,54	2121/rb	1,38	1000-1200	39.4 - 47.2
XL	2221/rb	1,60	2321/rb	1,66	2421/rb	1,66	2121/rb	1,50	1000-1200	39.4 - 47.2
S	2221/re	1,22	2321/re	1,25	2421/re	1,24	2121/re	1,10	1450-1700	57.1 - 66.9
M	2221/re	1,46	2321/re	1,51	2421/re	1,50	2121/re	1,36	1450-1700	57.1 - 66.9
L	2221/re	1,62	2321/re	1,65	2421/re	1,64	2121/re	1,48	1450-1700	57.1 - 66.9
XL	2221/re	1,72	2321/re	1,77	2421/re	1,76	2121/re	1,60	1450-1700	57.1 - 66.9
S	2221/f	1,16	2321/f	1,23	2421/f	1,22	2121/f	1,04	1200-1450	47.2 - 57.1
M	2221/f	1,48	2321/f	1,49	2421/f	1,48	2121/f	1,34	1200-1450	47.2 - 57.1
L	2221/f	1,64	2321/f	1,69	2421/f	1,68	2121/f	1,50	1200-1450	47.2 - 57.1
S	2221/bf	1,11	2321/bf	1,17	2421/bf	1,17	2121/bf	0,99	1000-1200	39.4 - 47.2
M	2221/bf	1,43	2321/bf	1,49	2421/bf	1,49	2121/bf	1,29	1000-1200	39.4 - 47.2
L	2221/bf	1,59	2321/bf	1,65	2421/bf	1,65	2121/bf	1,45	1000-1200	39.4 - 47.2
S	2221/ef	1,24	2321/ef	1,27	2421/ef	1,26	2121/ef	1,12	1450-1700	57.1 - 66.9
M	2221/ef	1,52	2321/ef	1,57	2421/ef	1,56	2121/ef	1,38	1450-1700	57.1 - 66.9
L	2221/ef	1,68	2321/ef	1,75	2421/ef	1,74	2121/ef	1,54	1450-1700	57.1 - 66.9
S	2221/rf	1,24	2321/rf	1,25	2421/rf	1,24	2121/rf	1,08	1200-1450	47.2 - 57.1
M	2221/rf	1,54	2321/rf	1,57	2421/rf	1,56	2121/rf	1,44	1200-1450	47.2 - 57.1
L	2221/rf	1,72	2321/rf	1,73	2421/rf	1,72	2121/rf	1,60	1200-1450	47.2 - 57.1
XL	2221/rf	1,86	2321/rf	2,03	2421/rf	2,02	2121/rf	1,72	1200-1450	47.2 - 57.1
S	2221/rbf	1,18	2321/rbf	1,24	2421/rbf	1,24	2121/rbf	1,02	1000-1200	39.4 - 47.2
M	2221/rbf	1,48	2321/rbf	1,54	2421/rbf	1,54	2121/rbf	1,38	1000-1200	39.4 - 47.2
L	2221/rbf	1,66	2321/rbf	1,72	2421/rbf	1,72	2121/rbf	1,54	1000-1200	39.4 - 47.2
XL	2221/rbf	1,80	2321/rbf	1,86	2421/rbf	1,86	2121/rbf	1,66	1000-1200	39.4 - 47.2
S	2221/ref	1,28	2321/ref	1,31	2421/ref	1,30	2121/ref	1,14	1450-1700	57.1 - 66.9
M	2221/ref	1,58	2321/ref	1,75	2421/ref	1,74	2121/ref	1,46	1450-1700	57.1 - 66.9
L	2221/ref	1,76	2321/ref	1,77	2421/ref	1,76	2121/ref	1,64	1450-1700	57.1 - 66.9
XL	2221/ref	1,90	2321/ref	1,93	2421/ref	1,92	2121/ref	1,78	1450-1700	57.1 - 66.9


\*Span width: Distance between the attachment points on the bed / size XXL on request

Tables: approximate measurements – about 8 % shrinkage possible (products are not pre-washed)



**SEGUFIX®-Magnetic locking system black**

Item	Consisting of:
1206	1 Magnetic Button black / segufix.com 1 Metal Pin
1208	1 Magnetic Key red
1204	1 Magnetic Button black / segufix.com
1205	1 Metal Pin
1209	1 Magnetic Key red with key ring




**SEGUFIX®-Patent Twist locking system green**

Item	Consisting of:
1306	1 Patent Twist Button green 1 Square Pin
1308	1 Patent Twist Key yellow
1304	1 Patent Twist Button green
1305	1 Square Pin

**Intended use of the Lock:**  
Closing of the SEGUFIX® systems.

**Intended use of the Key:**  
Quick release of the locks.

**Caution:** Each locking system only works with its own component parts and is not interchangeable with other locking systems.



**Push locking system**

Item	Consisting of:
1406	1 Push Button 1 Push Pin
1408	1 Push Lock Key
1404	1 Push Button
1405	1 Push Pin



**WARNING:** Read the Safety Information prior to using the SEGUFIX®-Locks and Keys.

- Do not submerge locks in any type of liquid. This can cause malfunction.
- Locks and keys must not be laundered.
- Test locks for proper function before each use.
- Examine the head of the pin for damage (wear and tear, rounding off or deformation).
- When closing the lock make sure that the button locks securely onto the pin. You should not be able to pull the button off of the pin.
- When opening locks, ensure movement between the upper and lower points of contact (item 1206). Therefore each lock can hold up to a maximum of 4 belt layers (exception: products with “r” or “f” in the item number have thicker belts. These allow for a maximum of 3 belt layers!). When using the Patent Twist Lock (item no. 1306) use no more then 4 belt layers. When using the Push Lock (item no. 1406) use no more then 3 belt layers.
- If damaged locks fail to open, straps can be cut with scissors or belt cutters.
- For more product info visit us online at [www.segufix.com](http://www.segufix.com) or [www.segufix-shop.com](http://www.segufix-shop.com)

## Instructions for use: SEGUFIX®-Locks and Keys

**Item 1206 / 1208** (SEGUFIX®-Magnetic locking system black)

**Item 1306 / 1308** (SEGUFIX®-Patent Twist locking system green)

**Item 1406 / 1408** (Push locking system)

### How to open the Locks:

#### Item 1206 / 1208

- Set the red Magnetic Key with a light push on top of the black Magnetic Lock.
- Lift both the Magnetic Key and the Magnetic Button off of the Metal Pin (fig. 1).



#### Item 1306 / 1308

- Set the yellow Magnetic Twist Key on top of the green Magnetic Twist Lock.
- Unlock the Patent Twist Lock by turning the Patent Twist Key, and then lift both the Patent Twist Key and the Patent Twist Button off of the Square Pin (fig. 3).



#### Item 1406 / 1408

- Push the two pins of the Push Lock Key into the holes on top of the Push Button to unlock.
- Lift both the Push Lock Key and the Push Button off of the Push Pin (fig. 2).



### How to close the Locks:

- Set the respective button with the opening onto the pin so that it locks in place.
- Now test the lock for a secure hold: You should not be able to pull the button off of the pin.

**Caution:** Each locking system only works with its own component parts and is not interchangeable with other locking systems.

**Follow SEGUFIX®-Safety Information and Warnings on pages 8 - 11**

The following instructions for use are only applicable to original SEGUFIX<sup>®</sup> products. Do not use the SEGUFIX<sup>®</sup>-Restraint and Positioning System until you have read and understood the instructions for use. Only use up-to-date instructions and keep them in a safe and readily available location. Product and instruction must match. Failure to use this product properly can cause serious injury or death. Your patient's safety depends on your using the SEGUFIX<sup>®</sup> system correctly.

## 1. Requirements for use

It is the user's responsibility to comply with all national laws. Always consider restraining protocols of your work place and current legislation. The instructions for use must be adhered to.

SEGUFIX<sup>®</sup> products may only be used by individuals with proper training and knowledge of the product. Legal provisions and doctor's orders must be strictly observed! The attending physician must ensure, that the patient's medical condition allows the application of the SEGUFIX<sup>®</sup> system. Watch for unlikely allergic reactions.

## 2. Intended use

SEGUFIX<sup>®</sup> products may only be used for the medical purposes, which they were designed for.

## 3. Remove dangerous objects

Remove all objects from the reach of the patient, which could lead to injury or may damage the belts (e.g. glasses, jewelry and sharp objects).

## 4. Inspect products before use



Inspect before each use for broken stitches, torn, worn out or frayed parts and damaged eyelets. Damaged products could malfunction and cause injury. Do not use damaged or worn out products.

## 5. When opening, locks may not be under tension

When opening locks, ensure movement between the upper and lower points of contact (item 1206). Therefore each lock can hold up to a maximum of 4 belt layers (exception: products with "r" or "f" in the item number have thicker belts. These allow for a maximum of 3 belt layers!). When using the Patent Twist Lock (item 1306) use no more than 4 belt layers. When using the Push Lock (item 1406) use no more than 3 belt layers.



## 6. Always attach Side Positioning Straps

The SEGUFIX®-Standard must always be used with the Side Positioning Straps attached. Side Positioning Straps prevent the patient from helicoptering in bed or from getting over the edge of the bed causing compression of the chest, possibly resulting in death.

In some cases it is possible that patients get over the edge of the bed even with the Side Positioning Straps attached and the side rails in the up position. In these cases additional restraint measures must be taken.

**Important:** Insert the Side Positioning Straps through the respective loops (according to instructions) before securing them to the bed belt.

## 7. Beds, chairs and stretchers


- Only use SEGUFIX® products on beds, stable chairs and stretchers that allow the product to be attached to, in accordance with the instructions.
- Products must be attached in such manner, that they cannot get out of place.
- Note that straps and moving mechanics of beds, stretchers or chairs (or alternating pressure pads and the like), may interfere with each other causing damage, loss of function and injury.
- The attachment points must be free of sharp edges.
- SEGUFIX® products do not prevent beds, chairs or stretchers from tipping over or from unintended movement.

Helpful hints:

European standards for beds: DIN EN 60601-2-52

European standards for wheelchairs: DIN EN 12182 and DIN EN 12183

## 8. All side rails must be in the up position

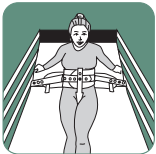

To prevent accidents always lock side rails in the up position (exception: 5-point-restraint).

**Caution:** Do not use with parted side rails except with a gap protector that effectively prevents the patient from slipping through the gap. An open gap could allow the patient to slip through even with the restraint fastened possibly causing compression of the chest.

## 9. Monitor patients in restraints



Patients in restraints require frequent monitoring.

Aggressive, agitated or restless patients and patients in danger of aspiration require constant monitoring.

## 10. Test closure-systems prior to each use



Check SEGUFIX®-Magnetic Locks black/segufix.com, SEGUFIX®-Patent Twist Locks green, Push Locks, Snap Locks, Velcro® Fasteners and Buckle Locks for a secure hold.

Do not submerge locks in any type of liquid. This may cause malfunction. If damaged locks fail to open, straps can be cut with scissors or belt cutters. Check the head of the pin for signs of wear and tear (rounding off

or deformation).

## 11. Pacemaker Warning



A distance of 10 cm (4 inches) must be maintained between a patient's pacemaker and magnetic locks and keys to prevent the pacemaker from possible malfunction (alternatively use Push Locks).

For possible interference of strong magnets with other patient implants please refer to manufacturer's instructions of these products.


## 12. Test products for correct attachment

Frequently check products for correct attachment. The straps must be tightly secured to the folding part of the bed, without folding part around the bed frame. The straps must stay in place tightly.

## 13. Only use original SEGUFIX® products

For safety reasons SEGUFIX® products must not be altered or be combined with third party products.

## 14. Use of SEGUFIX® products in vehicles

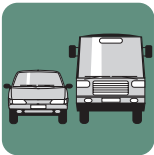

The SEGUFIX®-Transport Belt System does not replace the vehicle's legally required seat belt system. SEGUFIX®-Transport Belts may only be used as additional support or for restraint of a patient. Make sure that the patient can be rescued quickly in case of emergency (e.g. clearly visible and handy belt cutters).

## 15. Secure but not too tight

The restraints must be fastened tightly, but not interfere with breathing or blood circulation. Loosely fastened products may jeopardize the patient's safety. Frequently check products for correct position.

## 16. Storage

To prolong product life, store belts and locking systems in dry and dark conditions at about 20 °C.



## 17. Flammability of SEGUFIX®-Bed Restraints



According to DIN EN 71-2 section 4.3 non-flammable.

According to DIN 75200 the burning rate is zero.

However, products may not come in contact with cigarettes or be exposed to open flame, as this could permanently damage the products, and this could lead to injury.

## 18. Correct sizing

Always choose the right size is of essential importance. Make sure that the size matches body and waist sizes (refer to sizing tables). Products that are too small or too large compromise the patient's comfort and safety.

## 19. Disposal

Cardboard → recyclables

Plastic wrapping → recyclables or household garbage depending on local regulations

No longer functional products → household garbage

## 20. Declaration of conformity

For the declaration of conformity please visit:

[https://www.segufix.com/conformity\\_en.php](https://www.segufix.com/conformity_en.php)

## 21. Expert recommendation to prevent accidents

(Prof. Dr.-Ing. U. Boenick, TU-Berlin)

1. Always use Side Positioning Straps
2. Lock non-parted side rails in the up position  
(exception: five-point-restraint)
3. The use of additional restraint elements may be necessary to eliminate bed egress.

## 22. Ensure prior to each use

- Test locks for a secure hold
- Visual inspection of belts, straps and eyelets
- Ensure sufficient medical indication
- Check legal requirements

## 23. Serious incidents

Serious incidents must be reported to the responsible authorities and to the manufacturer.

# Care Instructions

- **SEGUFIX®-Restraints** (white) are made of a cotton/rayon mix and are washable in hot wash cycles up to 95 °C. Lower temperature and laundry nets will prolong product life. Avoid methods that might bend metal parts of the system (e.g. laundry line, dry pressing, rotary irons), eyelets may be damaged.  
Not pre-washed, up to 8 % shrinkage. Store in a dry place.  
Detergents: non-aggressive, no fabric softener, no bleach.  
Dryer: delicate cycle at low temperature and preferably in the laundry net.  
Disinfection: may be disinfected using chemothermal disinfection washing procedure. Overdosing can damage the product.
- **SEGUFIX®-Transport Belts** (coloured/black) are made of polyester and are washable in warm wash cycles up to 30 °C.  
Avoid methods that might bend metal parts of the system (e.g. laundry line, dry pressing, rotary irons), eyelets may be damaged. Store in a dry place.  
Detergents: non-aggressive, no fabric softener, no bleach.  
Dryer: delicate cycle at low temperature and preferably in the laundry net.
- **SEGUFIX®-Locks and metal pins:** remove all locks and metal pins before washing or drying. Laundering will damage locks and possibly cause them to fail. Locks may be cleaned and deinfected with a damp cloth. Never submerge locks in any type of liquid.
- **Velcro® Fastener:** fasten before washing and drying to prevent lint buildup. Occasionally use a stiff brush to remove lint from Velcro® Fastener. Even with normal use, Velcro® material has the tendency to collect lint. This can reduce the grip strength. Before each use, test grip strength of the Velcro® Fastener.
- **Durability:** depends on the intensity of use and is prolonged by gentle laundry procedures (detergent/laundry net).

---

Edition 19 · 2022-07

© by SEGUFIX®-Bandagen, 2022

The general business and delivery terms apply:

[https://www.segufix.com/PDFs/AGBs/Terms\\_and\\_Conditions.pdf](https://www.segufix.com/PDFs/AGBs/Terms_and_Conditions.pdf)

Product specifications subject to change. Deviations from pictures in colour and shape, errors and amendments excepted.

All rights reserved. Reprints or distribution using electronic systems/media, in its entirety or in parts, with prior explicit authorization in writing only.

All requests concerning this must be directed at:

## **SEGUFIX®-Bandagen**

Das Humane System GmbH & Co. KG

Allerbeeksring 33 · 21266 Jesteburg  
Germany

Phone: +49 (0) 41 83 500-0

Fax: +49 (0) 41 83 500-200

E-mail: [export@segufix.de](mailto:export@segufix.de)

[www.segufix.com](http://www.segufix.com)

[www.segufix-shop.com](http://www.segufix-shop.com)